


MDP

DOUBLE PAN BALANCE

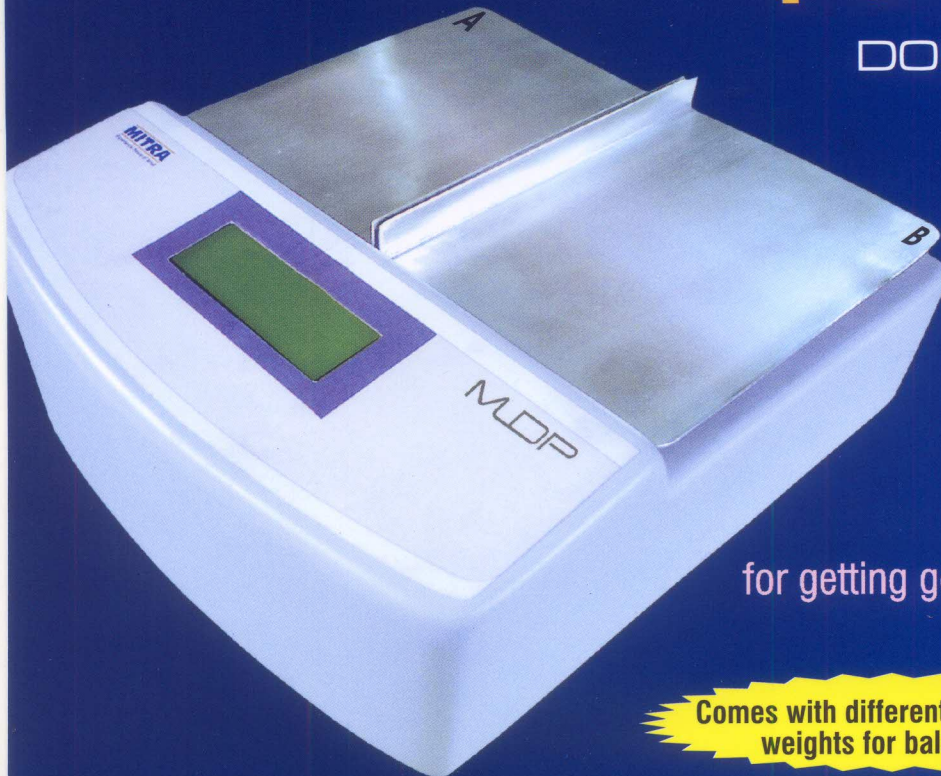

MDP is a micro controller based system which equally balance the blood bags after donation & before centrifugation for getting good quality of blood components.

Comes with different soft rubber weights for balancing

SPECIFICATIONS:

Load cell capacity	2500 grams
Weight sensor	Load cell
Weight display panel	Back light 4 line alphanumeric LCD display for better visibility indicates the weight difference of each pan
Indication Lamp	LCD display message after equal balance
Weight Difference display	LCD display with each pan weight difference
Calibration	Auto calibration
Accuracy (each pan)	Weight calibration tolerance : ± 2 gm.
Power source	220-240V AC 50Hz
weight	Approx 3 kg
Temperature characteristics	Operating : 0~30°C Storage : 0~40°C

For Further Information Please Contact:

Mitra Industries Pvt. Ltd.

14/4, Delhi-Mathura Road, Near Metro Station, Mewala Maharajpur, Faridabad, Haryana - 121 003. INDIA

Phone: 0129-2271068, 2271069. Toll Free No.: 1800 1020 307

e-mail: customersupport@mipl.in Web: www.mitraindustries.com